


Medical Instrument Cleaning Solution

This cleaning agent is specially designed for medical instruments, providing effective cleaning and disinfection. It removes contaminants and ensures the safety of medical equipment for patient care.

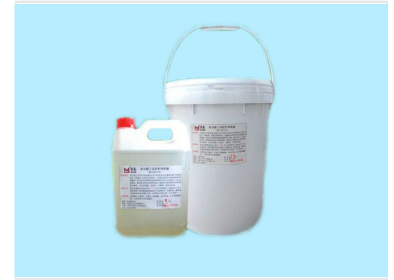

Product Overview

Professional Medical Instrument Cleaning Solution

This high-performance cleaning solution is specifically formulated to remove complex contaminants including oil, grease, and biological stains from medical instruments. Its super-concentrated, non-toxic formula ensures thorough sanitation without damaging equipment surfaces or leaving residue. Designed for versatility, it supports various cleaning methods such as ultrasonic, spray, and manual soaking to meet rigorous hygiene standards.

Technical Specifications

Optimal Cleaning Temperature

45 °C

Min Temp

75 °C

Max Temp

| | |
|------------------------------------|---|
| Compatible Cleaning Methods | Ultrasonic Cleaning, Spray Cleaning, Soaking, Manual Wipe |
| Recommended Concentration | 3% - 5% |

Safety & Compliance

Safety Profile

Heavy Metal Free • Phosphorus Free • Non-Toxic • Non-Corrosive

| | |
|---------------------------------|-----|
| Environmental Compliance | Yes |
|---------------------------------|-----|

Process Workflow

Standard Cleaning Workflow

- Prepare solution (3-5% concentration)
- Heat solution to 45°C-75°C
- Submerge parts into solution
- Apply ultrasonic, spray, or brush cleaning
- Rinse thoroughly with water
- Dry parts completely

Product Details

| | |
|-------------------|-----------|
| Shelf Life | 12 months |
| pH Value | 7 |